

# Introduction to Bulk RNAseq data analysis

## QC of raw reads with FastQC - Solutions

### Contents

1. Using FastQC to QC raw reads	1
Exercise 1 . . . . .	1

## 1. Using FastQC to QC raw reads

### Exercise 1

1. a) Check the location of the current directory using the command `pwd`
  - b) If the current directory is not `Course_Materials`, then navigate to the `Course_Materials` directory using the `cd` (change **d**irectory) command:

```
cd ~/Course_Materials
```

2. a) Use `ls` to list the contents of the directory. There should be directory called `fastq`
  - b) Use `ls` to list the contents of the `fastq` directory:

```
ls fastq
```

```
SRR7657883.sra_1.fastq.gz SRR7657883.subset_2M.sra_1.fastq.gz
SRR7657883.sra_2.fastq.gz Test_adapter_contamination.gq.gz.
SRR7657883.subset_2M.sra_2.fastq.gz
```

You should see two fastq files called `SRR7657883.sra_1.fastq.gz` and `SRR7657883.sra_1.fastq.gz`. These are the files for read 1 and read 2 of one of the samples we will be working with.

3. Run `fastqc` on one of the fastq files:

```
fastqc fastq/SRR7657883.sra_1.fastq.gz
```

This creates two files in the `fastq` directory. The first is the QC report in html format and the second is a zip file containing the data summary data used to generate the report. `> ⇒ SRR7657883.sra_1_fastqc.html`  
`> ⇒ SRR7657883.sra_1_fastqc.zip`

4. Open the html report in a browser and see if you can answer these questions:

## ✓ Basic Statistics


Measure	Value
Filename	SRR7657883_1.fastq.gz
File type	Conventional base calls
Encoding	Sanger / Illumina 1.9
Total Sequences	30120695
Total Bases	4.5 Gbp
Sequences flagged as poor quality	0
Sequence length	150
%GC	49

A) What is the read length?

**150**

B) Does the quality score vary through the read length?

## ✓ Per base sequence quality

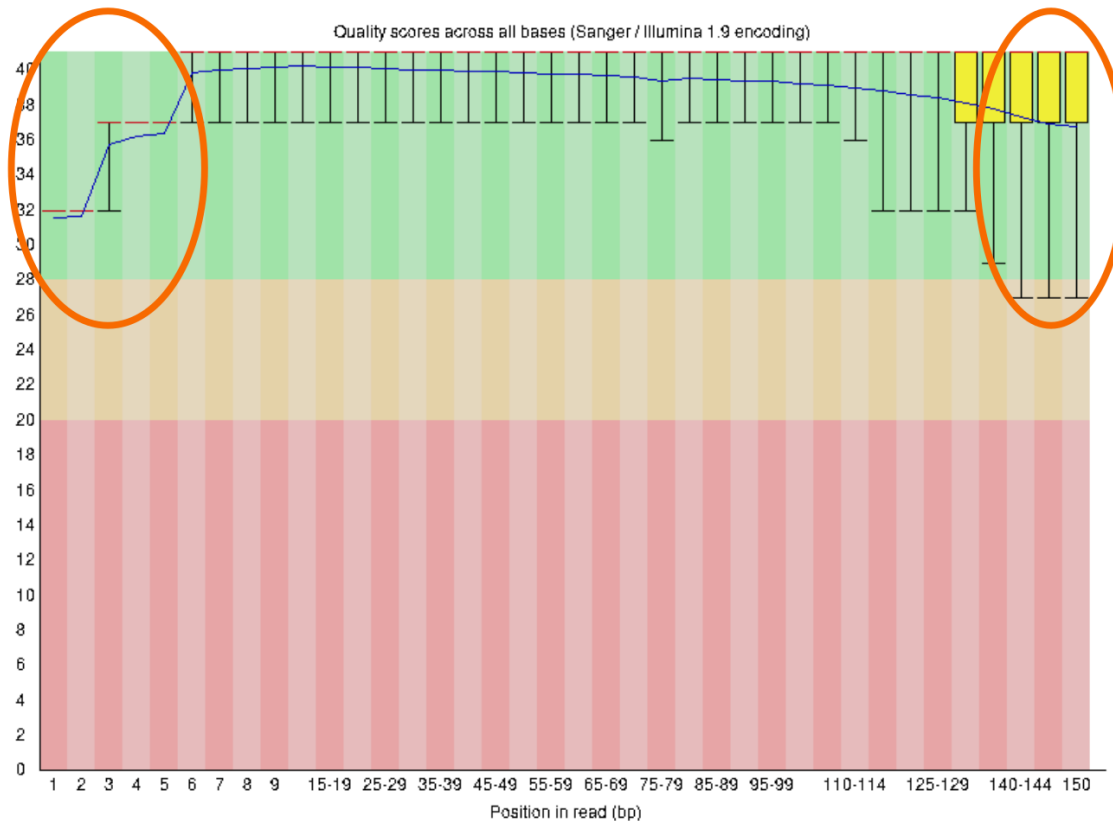

Yes, the first few bases and the last few bases are typically of lower quality.

C) How is the data's quality?

Overall, pretty good.